Presstech

MANUAL

ENGLISH	4
ESPAÑOL	32
FRANÇAIS	60
DEUTSCH	80
ITALIANO	116
PORTUGUÊS	144

PressTech

USER MANUAL

Congratulations on the purchase of your
SKINVITY PressTech device!

You have taken the leap into a new,
liberating, more enjoyable and efficient way
of personal care.

This manual describes the PressTech product,
how to use it and get the most out of it.

Additional information can be found at
[skinvity.com](https://www.skinvity.com) or by scanning this QR:

Index

Important Notice	8
Safety	8
Warnings	9
Precautions	10
What is PressTech	12
Accessories	14
Controls	15
Programs	16
Using PressTech	18
Recommendations for use	21
When using your PressTech	23
Duration of treatment	23
Frequency between treatments	24
Main benefits of PressTech	25
Troubleshooting	27
Customer service	28
Warranty	28

Important Notice

© SKINVITY
All rights reserved.

PressTech is not a medical device.

Design and technical specifications are subject to change without notice. Manufactured by FEMTECHNOLOGY SL.

All intellectual property rights in this publication are owned by FEMTECHNOLOGY SL. and are protected by applicable copyright laws and international treaty provisions.

FEMTECHNOLOGY SL. holds all rights not expressly granted. No part of this publication may be reproduced in any form or used to make any derivative work without the prior written approval of FEMTECHNOLOGY SL.

SAFETY

IMPORTANT: Be sure to read the entire user manual before using the product.

It is important to observe the WARNINGS and PRECAUTIONS when using PressTech.

WARNINGS

PressTech is an electrical device and it is important that you use it properly. Follow all safety WARNINGS and PRECAUTIONS:

- Do not touch the device if it falls into water. With dry hands, unplug it quickly.
- Do not place or store the device near any source of water, such as a sink or bathtub.
- Do not submerge the device underwater.
- If the device is plugged in, do not leave it unattended.
- Before each use, make sure that the system and its power supply show no visible damage. Do not use if you notice any damage to the power cord.
- Use the device only with the power cord and power supply provided.

- Keep the device out of the reach of children.
- Always use it with dry hands.
- After use or before cleaning, switch off the device and unplug it from the electrical outlet.
- Use the device only for its intended purpose as described in the manual.

PRECAUTIONS

The PressTech device is completely safe and generates no harmful side effects.

However, as a PRECAUTION, we recommend not to use it in the following cases:

- Children under 16 years of age.
- Persons with complex pathologies.
- In case of ulcers, wounds, bone fractures or risk areas.
- With fever or temperature above 37.5°C.

- If you suffer from cancer, heart disease or have a pacemaker.
- In case of hypertension, with osteosynthetic material or infectious diseases without medical supervision.
- In case of pregnancy, consult your doctor and never apply to the abdomen.
- In case of high blood pressure, we do not recommend mode 2.

At SKINVITY, we always recommend that you consult a specialist if you have any doubts.

What is PressTech

PressTech is a device used in aesthetic and self-care treatments to apply controlled pressure to different parts of the body. It aims to improve blood and lymphatic circulation, reduce fluid and toxin retention and alleviate issues such as tired and swollen legs.

It provides a comprehensive improvement in well-being.

Operating it is simple and intuitive, but it is essential to read the manual and know it well to get the most out of it.

Basic components:

1. **Air chamber (trousers):** This is the accessory that is placed around the limbs and reaches the abdomen and buttocks.
1. **Air compressor (control box):** This component generates the controlled air pressure that is used during the session.
1. **Tubing and connectors:** The tubing connects the air compressor to the boots, sleeves or girdles, allowing the air pressure to be distributed evenly.

Once the device is assembled and turned on, the air compressor begins to inflate the various chambers in a specific sequence.


Typically, pressure is applied in upward waves, from the feet towards the torso. The pressure and duration of each inflation and deflation cycle is controlled according to the mode selected.

Accessories

- | | |
|-------------------------------|--------------------------|
| 1. Control panel (compressor) | 5. Air chambers |
| 2. Y-connector | 6. Expansion accessories |
| 3. Hoses x2 | 7. Travel bag |
| 4. Power supply | |


Controls

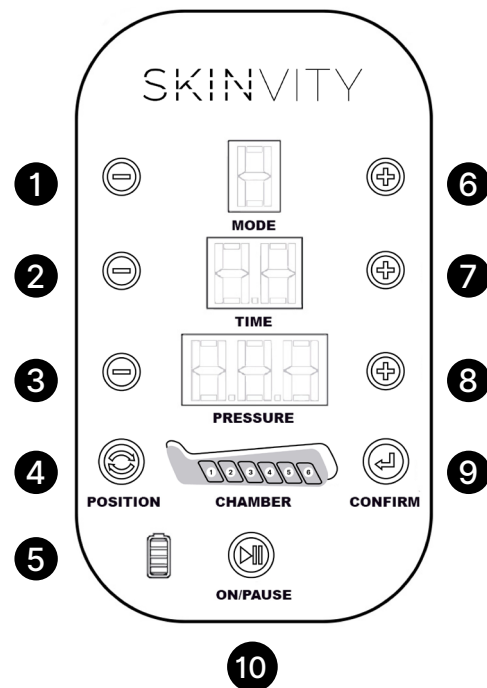- | | |
|----------------------|----------------------|
| 1. Previous mode | 6. Next mode |
| 2. Reduce time | 7. Increase time |
| 3. Reduce pressure | 8. Increase pressure |
| 4. Massage position | 9. Confirm position |
| 5. Battery indicator | 10. On/Pause |


Programmes

PressTech has 6 different compression modes that favour certain effects when used alone or in combination:

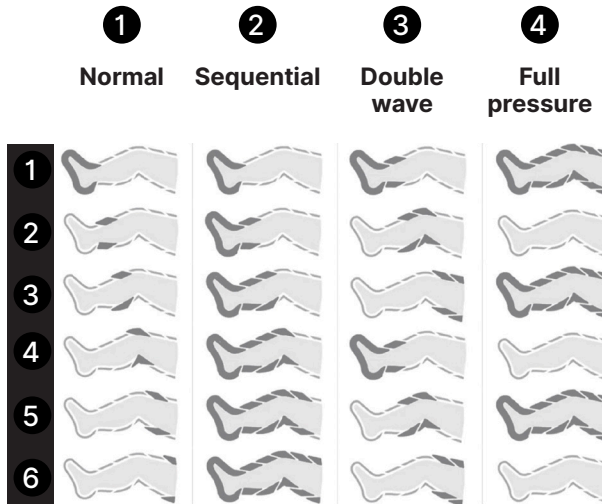- 1. Normal:** activates one chamber at a time, in order from the foot towards the torso. By inflating one and deflating the previous one, it offers an oxygenating and relaxing action for tired legs. Great draining capacity of liquids and toxins.

- 2. Sequential:** pressure exerted continuously on all the chambers, in sequence from the foot towards the torso, but maintaining the pressure. This treatment improves blood circulation.
- 3. Double wave:** the chambers press in groups of two to boost circulation. It also has a great draining action of liquids and toxins.
- 4. Full pressure:** all the chambers are activated and deactivated at the same time, contributing to a relaxing and effective action for drainage of fluid retention and the oxygenation of tissues and muscles.

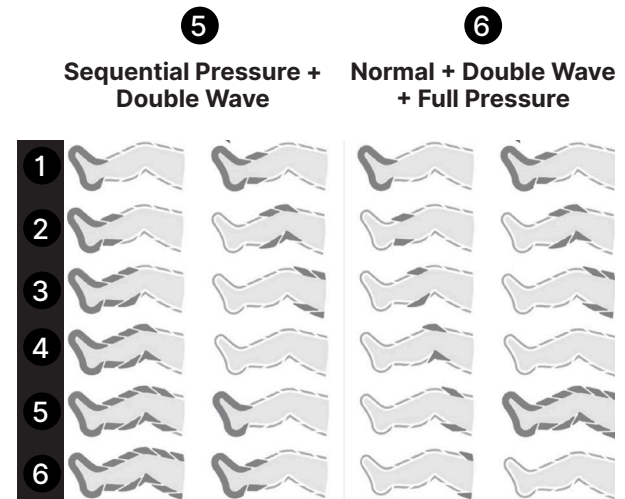5. Sequential Pressure + Double Wave: We alternate modes 1 and 2, so that initially only one is activated and then they are activated in groups of two. These two combined modes stand out for their great circulatory and draining action of liquids and toxins.

6. Normal + Double Wave + Full Pressure: Modes 1, 2 and 3 are combined, starting with the activation of a single chamber, then two by two, until all the chambers are activated at the same time.

Using PressTech

Before you get into your fabulous new PressTech device, follow these instructions:

1. Open the box, remove the device, air chamber, power supply, hoses and Y-connector.
2. Connect the power supply to the device and in turn to a power outlet.
3. Allow the device to fully charge.
4. Disconnect the power supply from the device.

When using your PressTech:

1. Connect the device to the air chamber (the trousers), using the hoses and the Y-connector, which should be at the top of the trousers, so that you can easily handle it. Check that there are no knots or bottlenecks.
2. If you need the extender accessory, which adds 5 cm on each side (size 42-44), attach it to the trousers. It's very easy, just add it using the side zips.
3. Put the trousers on in a comfortable place, either half lying or lying down, and make sure you pull them up well, to just below your chest.
4. Turn the device on using the side switch.
5. Select the duration of the programme using the "Time" buttons. It defaults to 30 minutes and can be easily raised or lowered from the controller.
6. Select the desired pressure using the "Pressure" buttons.

IMPORTANT: *Adjust this value until you get moderate pressure, it is best to start from less to add more.*

7. Select the desired programme (see pages 16 and 17) using the buttons in the “Mode” section. It is only possible to change the mode when the device is paused.
8. Use the “Position” button to select the chambers you wish to deactivate. Then confirm with the “Confirm” button. You cannot only massage a specific part of the legs or abdomen.
9. Press the “ON” button to start. Enjoy your self-care session from the comfort of your home.

To stop using it:

10. Press the “Pause” button.
11. Press the switch on the side of the device, leaving it in the 0 position.
12. Disconnect the hoses from the control unit; to avoid damage due to wear and tear, avoid disconnecting the hoses from the trousers. It is possible to store it in the bag with the hoses attached.

Disconnecting

- Avoid pulling on the hoses.
- To disconnect the hoses, do so directly from the connector.
- Press the push buttons on the side of the connector to facilitate disconnection.

Recommendations for use

Here are a few tips to help you get the most out of your brand-new pressure massager:

1. It is advisable to be hydrated before the treatment, but do not drink water or eat food just before starting.
2. To boost the elimination of localised fat, we recommend a manual massage beforehand (or with any instrument), to mobilise and warm up the fat. This will make its elimination with PressTech more efficient.
3. For a complete treatment with the aim of eliminating liquids, toxins and accumulated fat, we recommend starting with 10 minutes of mode 2 and then 10 minutes of mode 1, with steady pressure.

4. To finish with a relaxing treatment, add 5 or 10 minutes of mode 4 with very gentle pressure and/or 5 or 10 minutes of mode 3 with gentle pressure.
5. Before or after the treatment it is possible to apply cosmetics with firming, slimming or anti-cellulite effect.
6. It can be combined with other smart devices for home use.
7. We recommend not to exceed 30 minutes; no additional results are achieved.

When to use PressTech

PressTech is a device that uses air pressure to revitalise the body. It is apt for you if:

- You have cellulite or flaccidity.
- You are looking to improve inflammation and mitigate the feeling of swollen legs.
- You tend to retain fluids or toxins.
- Your legs ache.
- Presstech is not designed to perform lymphatic drainage or to reduce, eliminate or destroy adipose tissue under medical indication or prescription of health professionals.

Duration of treatment

The recommended usage time for our PressTech device may vary depending on the treatment objectives.

To help improve circulation and reduce fluid retention, shorter sessions may be more effective (15'- 20').

To improve the appearance of cellulite or to promote relaxation, longer sessions (30') are often more suitable.

Frequency between treatments

Generally, to improve cellulite and reduce fluid retention, the optimal frequency would be 1 to 2 times per week. Remember that you can combine it with other devices and cosmetics.

For maintenance, you should ideally use our PressTech device 2 to 4 times a month.

Main benefits of PressTech

- ✓ **Reduction of fluid retention and toxins:**
Our PressTech device uses sequential air pressure to apply compression and release to different parts of the body, which helps move fluid build-up in the tissues and improves lymphatic and blood flow. This can help reduce fluid retention and promote the elimination of toxins.
- ✓ **Improves blood circulation:** the sequential pressure applied to different parts of the body during a pressotherapy session helps to improve blood circulation, which is beneficial - and a relief - in combating the feeling of swollen, tired legs.
- ✓ **Cellulite reduction:** our PressTech device can help improve the appearance of cellulite by stimulating circulation and helping to reduce liquids and accumulated toxins, the presence of which is closely related to the appearance of cellulite on the skin.
- ✓ **Relaxation and well-being:** PressTech sessions are relaxing and can help reduce stress and tension. This can have a positive impact on overall wellbeing.

- ✓ **Improved skin appearance:** by improving circulation and reducing fluid retention, pressure therapy can contribute to more oxygenated, smoother and firmer skin.
- ✓ **Complement for weight loss and self-care:** although PressTech is not designed as a weight loss device, if you are following a care plan that includes diet and exercise it can be a useful addition and enhance its effects.

SKINVITY would like to remind you that it is important to note that the pressure massage technique can vary in terms of results and effectiveness from person to person, and is generally recommended as part of a holistic approach to body and skin care.

We must emphasise that too much pressure can be harmful, and can cause an uncomfortable feeling of tightness in the treated area, which can result in anxiety or discomfort during the session. It is necessary to go from less to more and stop when we feel that we are already receiving a firm massage, without becoming oppressive. You don't need to apply strong pressure to get the benefits of your PressTech device.

Troubleshooting

If your PressTech device is not working, follow these steps before contacting the customer service centre:

- Make sure the power cord is properly connected to the device and the battery is charged.
- Check that the side switch is in the "I" position and that the markers on the device are lit.
- Make sure that the ducts or hoses are not obstructed or choked.

Customer Service

If you have any problems with PressTech, please contact the retailer where you purchased it or visit skinvity.com, where you can read our FAQs.

Our ultimate aim is to provide femtech wellness to as many women as possible. To help us achieve this great objective and contribute to more women knowing and enjoying PressTech, you can send us your comments and suggestions, either through our social networks or by sending us an email to info@skinvity.com. We want you to take care of yourself.

Legal Warranty

For reasons of safety and hygiene, we do not accept returns of opened or used products, except for manufacturing defect.

In accordance with current regulations, all our products come with a three-year legal warranty from the date of delivery.

Warranty Conditions

- During the first two years, it is presumed that

any non-conformity detected already existed at the time of delivery unless SKINVITY proves otherwise. However, if necessary, SKINVITY may request additional information from the customer to assess the condition of the product.

- From the third year onwards, the consumer must demonstrate that the non-conformity already existed at the time of product delivery.

Warranty Exclusions

The warranty does not apply in the following cases:

- Damage caused by improper use, impacts, falls, or improper handling of the product.
- Natural wear and tear due to regular use of the product.
- Repairs or modifications not authorized by SKINVITY.
- Use in conditions not recommended in the user manual.

- Use of the product for commercial purposes instead of personal use
- Attempt to open, disassemble, or tamper with the device or its accessories

***Please take care of your PressTech device
and it will take care of you.***

**Name and address of the manufacturer:
Femtechnology S.L. C.I.F. B16846966**

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