
Silicone Led Mask

MANUAL

SKINVITY
tech

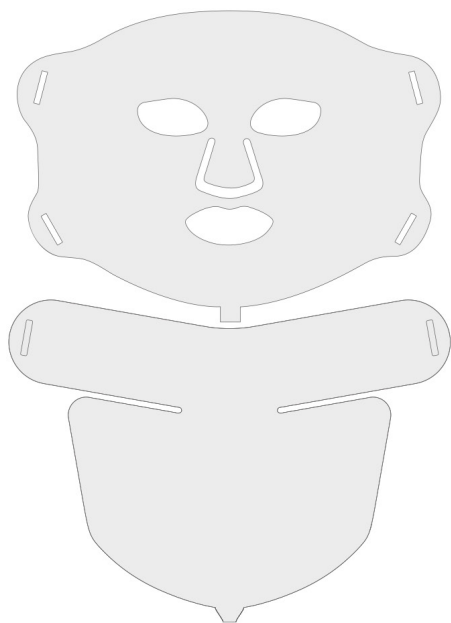ENGLISH	4
ESPAÑOL	28
FRANÇAIS	52
DEUTSCH	76
ITALIANO	100
PORTUGUÊS	124

Silicone Led Mask


MANUAL OF USER

Congratulations on the purchase of your
SKINVITY Silicone LED Mask device.

You have taken the leap into a new form of facial
care that is freer, more pleasant and more
efficient.

This manual describes the Silicone LED Mask
product and how to use it and get the most out
of it.

Additional information can be found on
skinvity.com or by scanning this QR code:

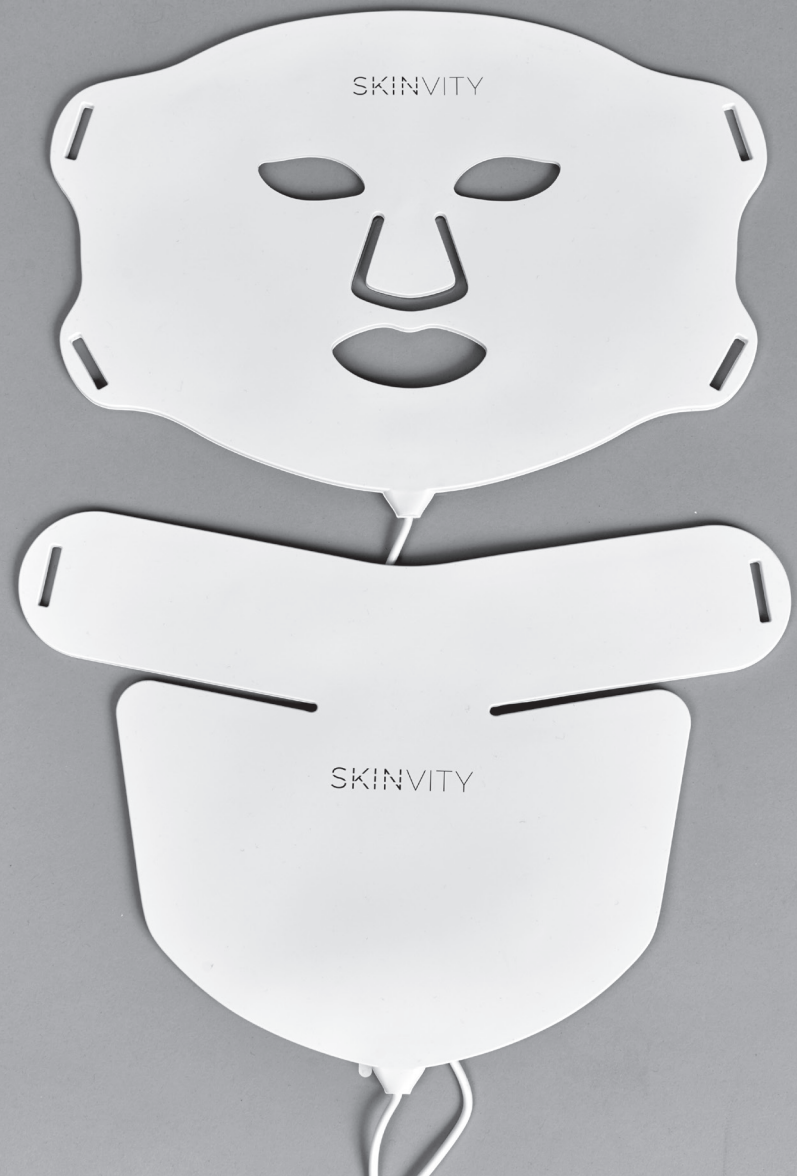Index

Important Notice	8
Safety	8
Warnings	9
Precautions	12
What is Silicone LED Mask?	13
Use of Silicone LED Mask	15
When to use Silicone LED Mask	18
Duration of treatment	18
Frequency between treatments	19
Main benefits	20
Troubleshooting	22
Customer service	22
Guarantee	23
Technical Specifications	25

Important Notice

© SKINVITY
All rights reserved.

Design and technical specifications are subject to change without notice. Manufactured by FEMTECHNOLOGY SL. All intellectual property rights in this publication are owned by FEMTECHNOLOGY SL. and are protected by applicable copyright laws and international treaty provisions.

FEMTECHNOLOGY SL. holds all rights not expressly granted. No part of this publication may be reproduced in any form or used to make any derivative work without the prior written approval of FEMTECHNOLOGY SL.

SAFETY

IMPORTANT: be sure to read the entire user manual carefully before using the product.

Always observe all WARNINGS and PRECAUTIONS when using your Silicone LED Mask.

Silicone LED Mask WARNINGS

Silicone LED Mask by SKINVITY is an electrical device that must be used correctly. To reduce the risk of injury from electric shock, observe all of the following WARNINGS and PRECAUTIONS:

- Do not touch the device if it falls into water. Unplug the device immediately with dry hands.
- Do not place or store the device near any source of water, such as a sink or bathtub.
- Do not immerse the device under water.
- When plugged in, do not leave the device unattended.
- Before use, make sure that the system and its power supply show no visible damage. Do not use if the power cord is damaged.
- Use the device only with the power cord and power supply provided.
- Keep the device out of the reach of children.
- Make sure your hands are completely dry when using the device.

- After your beauty routine or any maintenance, turn off the mask and unplug it from the electrical outlet.
- Use this device only for its intended purpose as described in this manual.
- The sensation during use of the device should be warm and pleasant. If you experience discomfort, discontinue use immediately.

Failure to follow the instructions in this manual could result in skin damage.

Unlike laser and IPL treatments, red and infrared and blue LED lights do not heat up the skin and do not cause pain, allowing the treatment to be performed without the need for a recovery period afterwards. These lights trigger reactions at the cellular level to naturally increase collagen production and promote tissue repair.

Charger WARNINGS:

- Do not disassemble or modify the battery.
- Do not expose to fire or extreme heat.
- Do not immerse in water or liquids.
- Do not pierce or crush the battery.
- Use only the original charger.
- Keep out of reach of children.

PRECAUTIONS


LED light has been clinically proven to be non-offensive and does not generate harmful side effects. However, as a precaution it is recommended not to use if:

- Children under 16 years of age.
- You currently have cancer or a history of cancer, especially skin cancer or premalignant moles.
- You have serious concurrent conditions.
- You are pregnant.
- You have any heat-stimulated disease within the treatment area or have an active skin condition within the treatment area, such as sores, eczema or a rash.
- You are taking photosensitising medication.
- You have an autoimmune disease such as lupus and porphyria or if you have epilepsy.

If you have any other health concerns, please consult your doctor before using Silicone LED Mask.


What is Silicone LED Mask?

Silicone LED Mask is a facial device used in aesthetic treatments. A mask made of soft, ergonomic silicone that uses red and blue LED lights and infrared to improve the firmness, health and appearance of the skin.

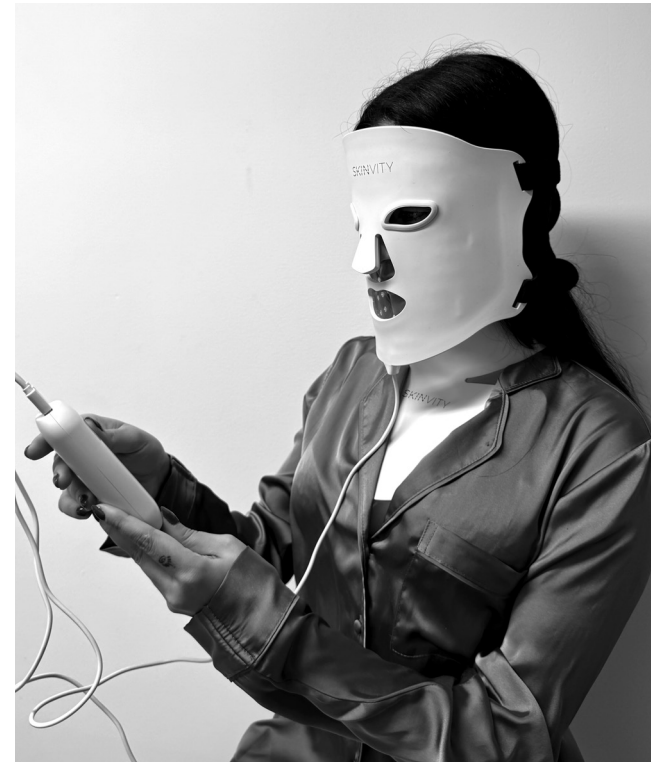
It is simple and intuitive to operate, but it is essential to get to know it well to get the most out of it:

Basic components:

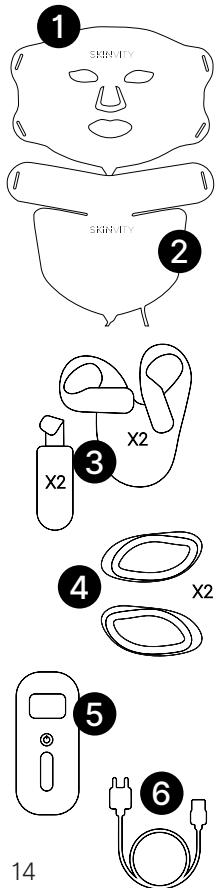1. Face accessory: the mask that is placed on the face.
2. Neck and décolleté accessory: an extra accessory to also treat the neck and décolleté.
3. Velcro: placed on the two face and neck attachments to hold them in place.
4. Eye shields: allow you to keep your eyes open during treatment. Includes a spare.
5. Remote control: charges and serves as a battery. Easy and intuitive to use.
6. Power supply and cable: connects the remote control to the mains to charge it.

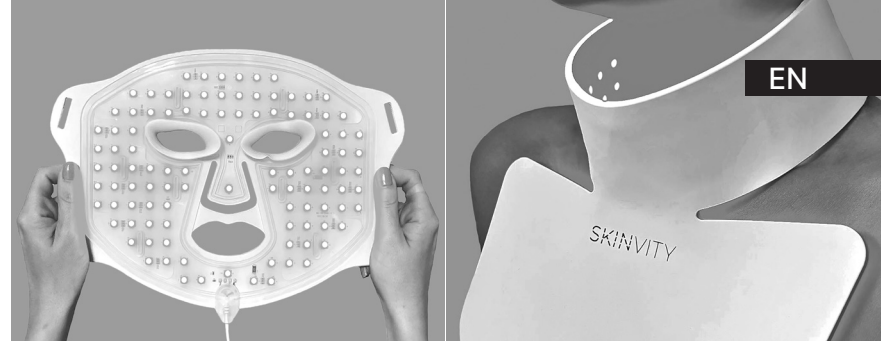Use of Silicone LED Mask

A few brief indications for a pleasant and efficient treatment:

1. Open the box containing the device. Remove the remote control, the device and its power supply.
2. Connect the power supply to the remote control of the device and, in turn, to a power socket.
3. Leave the device to charge until charging is complete.
4. Disconnect the power supply from the remote control.

We recommend using the mask with the battery charged, but not connected to the power supply. Using the Silicone LED Mask connected to the mains may cause the power to fluctuate during treatment.

When treating your skin with the Silicone LED Mask:

1. Make sure your skin is clean and dry.
2. Connect the controller to the mask.
3. Put on both elastic bands, comfortably adjusting the size 8.
4. Activate the remote control on the power button
5. Select the treatment mode according to the benefits you need using the "MODE" button.
6. Select the power level that will determine the time of your treatment on the "LEVEL" button.

The SKINVITY LED mask has 2 modes (red and infrared, and blue LED light) an 5 power levels, and the treatment time can range from 10 minutes for the highest power level to 30 minutes for the lowest power level.

Once you have chosen the desired treatment, relax close your eyes, and enjoy your moment of well-being.

To conclude:

1. Press and hold the power button for two seconds to turn the device off.
2. Disconnect the remote control from the mask and store both in the box.
3. Apply your usual facial treatment.
4. After using the mask, wipe gently with a cloth. Do not clean with chemical solvents. Do not get wet.
5. If the silicone tends to yellow after a while, wipe with an alcohol wipe.

When to use Silicone LED Mask

The SKINVITY Silicone LED Mask can be used at any time of the day, although we recommend using it in the evening, after cleansing the skin and before applying the products of the nightly facial routine.

Duration of treatment

Ideally, a shock treatment of 12 to 24 sessions (depending on the individual) should be carried out over a period of 6 to 8 weeks.

It is recommended to carry out 2 to 3 sessions every week until completion.

To maintain the effectiveness once the shock treatment has been completed, it is recommended to carry out 4 maintenance sessions per month.

The wonderful results on your skin will be maintained for around 9-12 months, although it depends on each person.

Frequency between treatments

Sessions may be performed more frequently than specified in this manual and the treatment may be applied for more than 30 minutes per day. Above 30 minutes each day at maximum power the treatment may not generate additional results.

If you opt for a 20 or 30 minute treatment session, our recommendation is to do two or three 10-minute treatments at power level 5, the optimum power level to ensure good results.

Main benefits of Silicone LED Mask

It is a revolutionary device, with visible results after just a few sessions.

The effects of treatment with red and infrared, and blue LED light translate into improvements in the tissues exposed to the light thanks to cell stimulation. Depending on the wavelength, we will have a colour with different applications, to comprehensively treat issues such as flaccidity, skin elasticity, increased dermis thickness, wrinkles and expression lines, pigmentation, sebum regulation, bacterial problems, marks, redness, inflammation and lack of luminosity, among others.

Following the treatment properly favors the production of collagen and elastin, improves skin tone, oxygenates the tissues and favors the absorption of cosmetic products. You can obtain a cared for and radiant skin.

Efficacy tests confirm the anti-wrinkle and firming efficacy of the Silicone LED Mask. After a study of 56 days of continuous use with 20 users, a test made by an independent laboratory confirms that:

- ✓ There is a significant increase in firmness of 30,8%.
- ✓ There is a significant increase in elasticity of 19%.
- ✓ There is a decrease in discolouration of 18,6%
- ✓ There is an increase in dermal thickness of 12%.
- ✓ There is a decrease in wrinkles in the nasolabial fold of 20,3%, in the forehead of 15,6%, between the eyebrows of 12% and 5,5% in the eye contour.

It is a non-invasive technique, an alternative to invasive aesthetic skin treatments.

It does not use light or ultraviolet (UV) lasers. Only LED and infrared light.

No side effects, does not cause inflammation. Painless and relaxing technique. No need for recovery.

Troubleshooting

If your Silicone LED Mask device is not working, please follow the following steps before contacting technical support:

- Check that the power cord is securely connected to the device and the battery is charged.
- Press the power button for a few seconds and verify that the remote turns on.

Customer service

If you have any questions or problems with Silicone LED Mask please visit skinvity.com, where you can check the product's frequently asked questions.

Our ultimate aspiration is to bring femtech wellness to as many women as possible. To help us achieve this great goal and help more women know and enjoy Silicone LED Mask, you can send us your comments and suggestions, both through our social networks and by sending us an email to info@skinvity.com. We want you to take care of yourself.

Legal Warranty

For reasons of safety and hygiene, we do not accept returns of opened or used products, except for manufacturing defect.

In accordance with current regulations, all our products come with a three-year legal warranty from the date of delivery.

Warranty Conditions

- During the first two years, it is presumed that any non-conformity detected already existed at the time of delivery unless SKINVITY proves otherwise. During this period, the consumer does not need to provide proof of the defect.
- From the third year onwards, the consumer must demonstrate that the non-conformity already existed at the time of product delivery.

Warranty Exclusions

The warranty does not apply in the following cases:

- Damage caused by improper use, impacts,

falls, or improper handling of the product.

- Natural wear and tear due to regular use of the product.
- Repairs or modifications not authorized by SKINVITY.
- Use in conditions not recommended in the user manual.
- Use of the product for commercial purposes instead of personal use
- Attempt to open, disassemble, or tamper with the device or its accessories

Please, take good care of your Silicone LED Mask Silicone LED Mask and it will take care of you.

Name and address of the manufacturer:
Femtechnology S.L. C.I.F. B16846966

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Technical Specifications:

Power supply

Model:	HDY 18650
Battery Type:	Li-ion Recargable
Rated voltage:	3.7V DC
Capacity:	2000mAh
Power:	7.4Wh
Maximum charging voltage:	4.2V DC
Recommended charging current:	1000mA (1A)
Maximum charge current:	2000mA (2A)
Final discharge voltage:	3.0V DC
Charging temperature:	0°C to 45°C
Discharge Temperature:	-20 °C to 60°C
Storage temperature:	-20°C to 40°C (up to 3 months), -20°C to 25°C (up to 1 year)
Dimensions:	18.5mm (diameter) x 65.5mm (height)
Weight:	Approx. 43.8g

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