Infrared Led Mat

MANUAL

SKINVITY

ENGLISH

4

ESPAÑOL

24

Infrared Led Mat

SKINVITY USER MANUAL

Congratulations on the purchase of your SKINVITY Infrared LED Mat device.

You have taken the leap into a new form of body care that is freer, more pleasant and more efficient.

This manual describes the Infrared LED Mat product and how to use it and get the most out of it.

Additional information can be found on skinvity.com or by scanning this QR code:

Index

Important Notice	6
Safety	6
Infrared LED Mat Safety Warnings	7
Precautions	10
What is Infrared LED Mat?	11
Use of Infrared LED Mat	13
When to use Infrared LED Mat	16
Frequency and duration of treatment	16
Main benefits	17
Troubleshooting	19
Customer service	20
Legal Warranty	20
Warranty conditions	21
Warranty exclusions	21
Technical Specifications	28

Important Notice

© SKINVITY

All rights reserved.

Infrared LED Mat is not a medical device.

Design and technical specifications are subject to change without notice. Manufactured by Femtechnology S.L. All intellectual property rights in this publication are owned by Femtechnology S.L. and are protected by applicable copyright laws and international treaty provisions.

Femtechnology S.L. holds all rights not expressly granted. No part of this publication may be reproduced in any form or used to make any derivative work without the prior written approval of Femtechnology S.L.

SAFETY

IMPORTANT: be sure to read the entire user manual carefully before using the product. They will help you understand how the product works and achieve the best possible results.

Always observe all WARNINGS and PRECAUTIONS when using your Infrared LED Mat.

Infrared LED Mat SAFETY WARNINGS

Infrared LED Mat by SKINVITY is an electrical device and it is important to use it properly. Follow all safety WARNINGS and PRECAUTIONS:

- Keep the device out of the reach of children.
- Do not leave the device unattended while it is plugged in.
- Never place the device where it might fall or be pulled into water or other liquids. Do not use it in high-humidity areas, such as the bathroom or shower.
- Do not use the mat if your skin is red or swollen from a previous treatment, or if it is in pain or sunburnt.
- Do not use the device if you have had a recent surgery or other surgical procedures.
- No modification of this device is allowed.

- Follow the recommended treatment time. Do not use the device excessively (overtreatment).
- Do not use the device if it is damaged or if its accessories (such as the cable or adapter) appear damaged. In that case, return it to SKINVITY for inspection and repair.
- Using different accessories from those specified or provided, could increase electromagnetic emissions or decrease electromagnetic immunity.
- Keep the equipment away from devices that emit frequencies such as control base radio as this may cause interference and unexpected treatment results.
- There is a potential risk of electromagnetic or other interference between this device and other equipment. If you suspect interference, move the Infrared LED Mat away from sensitive devices.
- Inspect the device, and do not use the device if you see any signs of internal battery leaks. In that case, contact our SKINVITY customer service.

- The sensation during use of the device should be warm and pleasant. If you experience discomfort, discontinue its use immediately.
- Use the device only with the power supply provided and supplied adapter.
- Do not attempt to open or repair the device or power supply.
- After your treatment, switch off the device and unplug it from the electrical outlet.
- Use the device only for its intended purpose as described in the manual.

Failure to follow the instructions in this manual could result in skin damage.

If you experience any adverse reactions to this device, stop using it immediately and consult your doctor.

If you feel any discomfort, such as skin discoloration, skin damage, or if any other concern arises, stop using the device immediately and consult your doctor.

PRECAUTIONS

Infrared and red LED light has been clinically proven to be non-offensive and does not generate harmful side effects. However, as a precaution it is recommended not to use it if:

- Children under 16 years of age.
- If you have a Type V or Type VI skin tone on the Fitzpatrick scale.
- If you have photosensitive disorders or genetic eye conditions.
- If you are pregnant.
- Do not use the device on infants or children
- Serious medical conditions, such as photosensitive (light-triggered) seizures, cancer or cancerous lesions, or having a pacemaker.
- If you are taking photosensitive medications.
- Any heat-stimulated disease in the treatment area or active skin condition (sores, eczema, rash). In

this case, ask your doctor before using Infrared LED Mat.

- If you have autoimmune diseases such as lupus or porphyria.
- If you have recently undergone surgery or other surgical procedures.

If you have any other health concerns, please consult your doctor before using Infrared LED Mat.

What is Infrared LED Mat?

The Infrared LED Mat is a body device used in aesthetic and self-care treatments, that emits 630nm, 850nm, 930nm and 1072nm lights, using four types of LEDs to reach this function.

The Mat is made of soft neoprene and is designed to safely warm the skin and underlying tissues through red and infrared LED light.

It aims to enhance skin revitalisation, refine texture and tone, and promote a soothing sense of relaxation and comfort, contributing to an overall


feeling of wellbeing. In addition, it features an extra pulsed magnetic field function to further enhance its benefits.

The device is simple and intuitive to operate, yet it is essential to become familiar with it to maximise its effectiveness

Basic components:

- **The Mat:** device designed for the user to lie inside covering the back and lower body.
- **Adapter:** providing the correct electrical power to the device. Use only the supplied adapter to avoid damage or automatic shutdown.
- **Storage bag:** to safely store and transport the device and its accessories, protecting them from dust, dirt, and damage.

Once the device is assembled and turned on, begins to emit heat, creating a pleasant sensation, it is not supposed to be uncomfortable.

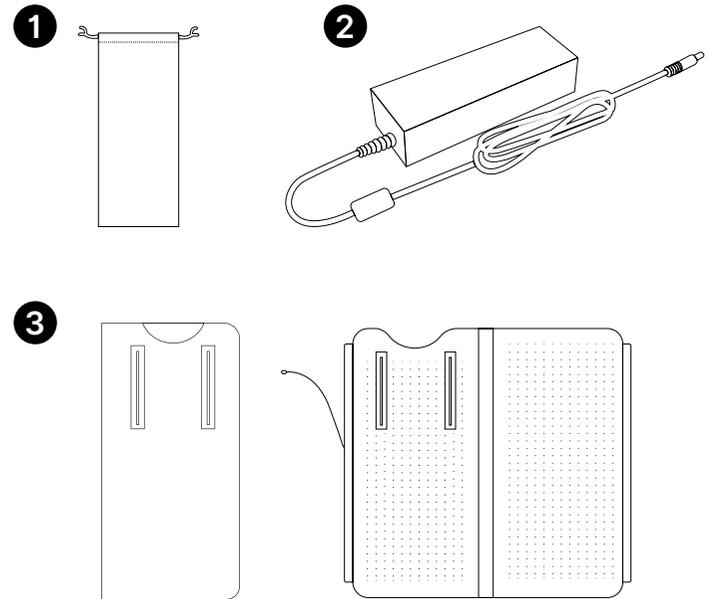1. Storage bag 2. Adapter 3. Infrared LED Mat

Use of Infrared LED Mat

A few brief indications for a pleasant and efficient treatment:

- Open the box containing the device. Remove the Mat, cable and its adapter.

- Ensure that the Mat, cable and adapter show no visible damage.
- Connect the adapter to the Mat.
- Connect the adapter to a power outlet.
- Place the Mat on a flat, comfortable surface where you can lie down.
- The power indicator light will begin to glow when the device turned on.
- Disconnect the power supply from the adapter.

The Mat must be connected to its adapter and, consequently, to the power supply. It is mandatory to use the supplied adapter.

When treating your skin with the Infrared LED Mat

- Make sure your skin is clean and dry.
- Clean for the device prior to use, unplug it from any power source.
- Ensure the Mat is completely dry before use.

- Place the Mat on a flat, dry surface before starting the treatment.
- Lie on the Mat, press and hold the power button for 2 seconds to start.
- Do not secure too tightly.

To conclude:

- Turn off the device by pressing the power button for 2 seconds, or will automatically shut down after 20 minutes.
- Disconnect the adapter and cable from the Mat and keep in storage bag after use.
- After its use, wait for the device to cool, do not over-treat.
- Clean the Mat gently with a damp cloth. Do not clean with chemical solvents.
- Do not get wet.
- Ensure the Mat is completely dry before storing in its bag after use.
- Store in environment with a temperature of -20~60°C and humidity of 30-90%RH.

For the first use, we recommend to test the device on one hand only to check light- sensitiveness. If you experience pain, discomfort or redness lasting more than 2 hours, please contact your doctor for advice.

When to use Infrared LED Mat

The SKINVITY Infrared LED Mat can be used at any time of the day, although we recommend using it in the morning or evening, after cleansing the skin and before applying the skincare products.

Frequency and duration of treatment

To achieve optimal results, we recommend using the device at least 2 to 3 times per week, although it can be safely used daily. The suggested treatment duration is 20 minutes per session.

Sessions longer than 30 minutes (approximately 1.5 treatments) do not provide additional benefits, though they are not harmful.

For maintenance, the device may be used 1 to 2 times per week.

Main benefits

An innovative device designed to support relaxation, comfort, ease of movement, and an overall sense of wellbeing by gently emitting red and infrared LED light over the body. Its main benefits are:

1. Improvement in skin firmness, elasticity and overall tissue quality: Red and infrared LED light stimulate cellular activity in the exposed areas, helping improve the appearance of body skin over time. This gentle stimulation supports natural collagen and elastin processes, contributing to a firmer, more elastic look, a smoother surface and an overall improvement in tissue density and suppleness across larger body areas.
2. Visible enhancement of skin appearance, tone and luminosity: Consistent use of the Infrared LED Mat can help refine skin texture and promote a more even-looking tone, improving the visible appearance of marks, uneven pigmentation and lack of radiance. By supporting oxygenation and revitalisation of the skin, the treatment helps achieve a brighter, more uniform and healthier-looking appearance, while also enhancing the absorption of cosmetic products applied afterwards.

3. Easing sensations of stiffness: Helps reduce everyday sensations of stiffness, promoting a greater feeling of flexibility and ease of movement.
4. Muscle relaxation: Encourages the body to feel more relaxed and less tense, supporting overall comfort during use.
5. Support for local circulation: May temporarily assist with circulation in the areas in contact with the device, contributing to a sense of local wellbeing.
6. Soothing effect on the skin: The gentle emission of red and infrared LED light provides a comforting and calming sensation on the skin.
7. Enhanced sense of wellbeing: Frequent use of the Infrared LED Mat can promote an overall feeling of comfort and relaxation, making it a valuable addition to body-care and wellness routines.

SKINVITY would like to remind you that it is important to note that in terms of results and effectiveness may vary from person to person, and is generally recommended as part of a holistic approach to body and skin care.

Troubleshooting

Trouble	Reason for troubleshooting	Solution
The device is not turning on	Ensure that the device is fully connected to the adapter	<ul style="list-style-type: none"> • If yes, you will see the power indicator light is on and solid • If not, please change the device
	Is the power indicator light illuminated after connecting?	<ul style="list-style-type: none"> • If yes, use it usually • If not, it may be damage
	The device is damaged	Please get in touch with info@skinvity.com
The device is not charging	Is the power cord firmly inserted into the device?	<ul style="list-style-type: none"> • Reconnect the charging cable and try again • Try a different wall outlet to determine if there is power to the outlet. • Please get in touch with info@skinvity.com
The device shuts off abruptly before the completion of the 20-minute treatment:	It may shut down because it overheated due to repeated use	<ul style="list-style-type: none"> • Wait for the device to cool, do not over-treat • After the device cools down, turn it on again. If the situation still occurs, please, get in touch with info@skinvity.com
The device shuts off abruptly after contacting highvoltage	This is the high voltage protection of the device itself, and it can be used after restarting	If the device cannot turn on, please, get in touch with info@skinvity.com
The device cannot be charged	Ensure that you are using a right adapter	If not, the device will not be able to charge

Customer service

If you have any questions or problems with Infrared LED Mat please visit skinvity.com, where you can check the product's FAQs.

Our ultimate aspiration is to bring femtech wellness to as many women as possible. To help us achieve this great goal and help more women know and enjoy Infrared LED Mat, you can send us your comments and suggestions, both through our social networks and by sending us an email to info@skinvity.com.

We want you to take care of yourself.

Legal Warranty

For reasons of safety and hygiene, we do not accept returns of opened or used products, except for manufacturing defect.

This product is covered by the statutory warranty of conformity in force in the country of purchase.

In case of malfunction, the device must be inspected exclusively by an authorized and certified service technician from Femtechnology S.L.

Femtechnology S.L. will not cover damages caused by conditions of use, storage and maintenance of the customer. Please take good care of your Infrared LED Mat.

Warranty conditions

This product includes a 3-year warranty covering manufacturing defects. The warranty is valid only with proof of purchase and does not cover damage caused by misuse, accidents, improper handling, or unauthorized modification.

Warranty exclusions

The warranty does not apply in the following cases:

- Damage caused by improper use, impacts, falls or improper handling of the product.
- Natural wear and tear due to regular use of the product.
- Repairs or modifications of the product not authorized by SKINVITY.

- Use in conditions not recommended in the user manual.
- Use of the product for commercial purposes instead of personal use.
- Attempt to open, disassemble, or tamper with the device or its accessories.

Please, take good care of your Infrared LED Mat and it will take care of you.

**Name and address of the manufacturer:
Femtechnology S.L. C.I.F. B16846966**

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